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### II. USE

Follow the instructions of installation in chapter indicated II maintenance.

The equipment has been planned to the aim to sterilize the suitable medical devices by means of warm air.

The equipment isn't usable in contact of the patient. The use is strictly recommended by doctors or nurses. It's necessary to keep the equipment where only authorized people can go.

The use of for any other applications, like the heating of alimentary substances or substances chemistries not in the specific watertight containers, strongly is advised against: the constructor cannot in some way to be thought in charge of the consequences of a not suitable use and not specifically previewed in this handbook.

**LOAD:** The trays are designed to withstand a maximum load of 25 kg loads in excess may affect the stability of the tray and the equipment, creating a danger of falling objects which could be at a high temperature.

**WARNING:** The items in the sterilizer and the inner walls of that can reach very high temperatures. It is recommended for safety and effectiveness of the sterilization cycle, to lock the door of the sterilizer and remove the keys when it is running.

In case of accidental contact with hot surfaces immediately contact a doctor.

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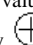
extension a value associated to the symbol indicated from the led illuminated to its left: the symbol  indicates the time and  indicates the temperature. Pushing key  change the visualization between time and temperature.

The programming of a cycle happens simply modifying the values of time and temperature second the application and the indications of the pharmacopoeia. This operation can be executed single on the value visualized to display. In the first instance, pushing

the key  , the value is chosen to modify (  time or  temperature); subsequently the value can be modified pressing the key 

in order to increase it or the key  in order to decrease it. Keeping pressed the keys, the variation happens in rapidly. Releasing the modification keys, the displayed value blink: the new value is saved when the blink stops (3÷4 second).

The time can be set up between 0 and 240 minutes, or to continuous cycle: keeping pressed the key  for some second after the attainment of the maximum value (240 minutes), the display show "inf" and

the led  blink; in this way the sterilizer maintains the temperature until when it turned off or reset from the customer. The temperature can be set up from 120°C to 200°C. Concluded the operation of modification of the parameters of the program (cycle), the equipment is ready for the use.

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Push the key to start the cycle. In the earlier stage display the extension "hea" written one: it indicates that the temperature has not still caught up 120°C (the minimal value that can be showed). The time set up for the cycle begins to diminish only when the equipment has caught up the set up temperature.

The Led  indicates that the resistance is in function. The cycle can be interrupted turning off the equipment or pressing the keys

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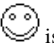
at the same time (reset). Interrupting the cycle, the display it visualizes an error message, to indicate that the contained material inside is not sterile (the message can be removed only executing a reset). This consideration is not worth for the infinite cycle. The classic configuration of dry sterilization cycle is 180°C for 30 min. The cycle is valid (material sterile) if display the extension the written “END” and the green led  is on.

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**ATTENTION:** If the cycle stop before “END” it's necessary to reset the cycle and restart.

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#### II.c. ERROR MESSAGES

During the cycle they can be verified of the marked anomalies from specific messages of error and the

lighting of the led red .

Tlo The temperature has not caught up the value setting up in a reasonable time: the probe or the resistance can be broken off.

Thi The temperature has exceeded excessive the set up value: it can be broken off the probe or the card of control.

Ert The temperature is come down quickly: Generally the door has been opened during the cycle.

OPD It marks them of participation of the security thermostat: probable breakdown of the electronic control.

ERR Cycle not finished: probably black-out an electrical worker or one manual interruption of the cycle.

Cancel the error messages by pressing the keys “+” and “-” simultaneously.

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### III. MAINTENANCE III.a. WARNING

**ATTENTION:** If the equipment was used in a manner not specified by the manufacturer, the protection provided may be compromised.

**WARNING:** The items in the sterilizer and the inner walls of that can reach very high temperatures. It is recommended for safety and effectiveness of the sterilization cycle, to lock the door of the sterilizer and remove the keys when it is running. In case of accidental contact with hot surfaces immediately contact a doctor.

**ATTENTION:** Do not place the unit on surfaces that can release smoke or flames when hot objects falling from the apparatus.

**ATTENTION:** Before introduce any object in the chamber be sure that the previously cycle is ended et the temperature of the inside walls is not to a dangerous level.

**ATTENTION:** When the cycle ended, before remove any object inside verify that his temperature isn't to a dangerous level.

**ATTENTION:** Do not open the door during the cycle. Keep out the key.

**ATTENTION:** If the cycle stop before "END" it's necessary to reset the cycle and restart.

**ATTENTION:** Before cleaning and maintenance, disconnect the electrical worker and wait the cooling of the metallic ones.

**ATTENTION:** In case of destruction follow the current requirements for the disposal of waste.

**ATTENTION:** The material subordinate to a cycle can be considered sterile single that you respect the parameters of the Official Pharmacopeia.

**ATTENTION:** Not to insert any thermo labile devices if not to the previewed maximum temperatures for sterilization of same and the always second Official Pharmacopeia.

**ATTENTION:** the constructor cannot in some way to be thought in charge of the consequences of a not suitable use and not specifically previewed in this handbook.

**ATTENTION:** The manufacturer can not in any way be held liable for any damage to the machine or other machines or persons resulting from an intervention that has not been executed by an authorized technician.

### III.b. INSTALLATION

Extract the equipment from the pack carefully. Keep the pack.

The equipment must be installed on one stable flat, far from thermo labile equipments and objects (not resistant to advanced temperatures to

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50°C) and in proximity of one taken of net electrical worker to 230V 50/60Hz supply of grounding. The equipment absorbs from 650 W to 1100W (depend of the model): it is advised to control the system electrical worker before connecting it and to use one taken exclusively for this application.

Maintain an adequate space for ventilation, leaving at least 10 cm between the casing of the machine and any adjacent walls or objects.

The equipment cannot be dipped in liquids neither bathe with liquid substances.

It is important that it comes at a distance held from inflammable or explosive substances.

**ATTENTION:** Do not place the unit on surfaces that can release smoke or flames when hot objects falling from the apparatus

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Periodically clean up the inside of the chamber with one sponge bathed with alcohol, do not use water spray do not throw directly water inside or outside of the machine.

Periodically control the conditions of the feeding cable and replacing it in case of cuts are observed.

**REPLACEMENT OF FUSES:** Make sure first that the device is turned off and unplugged from the electricity supply. Then opening the appropriate door located next to the socket on the left side of the machine, pull the fuses and replace them with suitable ones listed in Table 1 page. 4th

**REPAIR:** If the equipment showed abnormalities not directly solvable by following the instructions, you should contact an authorized dealer. The manufacturer can not in any way be held liable for any damage to the machine or other machines or persons resulting from an intervention that has not been executed by an authorized technician.

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## VII. GUARANTEE

### **MANUFACTURER'S LIMITED GUARANTEE**

LIMITED GUARANTEE: “M.B. PRO.MEC. SRL” O (M.B.) GUARANTEES THAT THE DRYSTER ST 20/40/60/80 LT IS FREE FROM DEFECTS OF MATERIAL AND PROCESSING FOR ONE YEAR FROM THE PURCHASING DATE.

THE GUARANTEE IS VALID ONLY IF THE PRODUCT WAS BOUGHT SOUND FROM M.B. OR FROM AN AUTHORIZED RETAILERS.

THIS GUARANTEE ISN'T VALID IF THE PRODUCT WAS USED IN AN ILLICIT WAY, IF WAS ALTERED, CHANGED OR REPAIRED BY PERSONNEL NOT AUTHORIZED BY M.B. OR BY UNAUTHORIZED SERVICE AGENTS.

IN THE CASE THAT THE PRODUCT SHOWS, IN ONE YEAR FROM THE PURCHASING, DEFECTS OF MATERIAL OR PROCESSING, IT MUST BE SENT TO M.B. OR TO A SERVICE CENTER AUTHORIZED BY M.B., WITH TRANSPORT EXPENSES PAID, ENCLOSING NAME, ADDRESS, THE PURCHASING DATE'S PROOF AND A SHORT DESCRIPTION OF THE NOTICED DEFECT.


M.B. AT ITS DISCRETION WILL MAKE THE REPAIR OR THE SUBSTITUTION OF THE FAULTY PRODUCT WITHOUT ANY EXPENSE BY THE PURCHASER.

POSSIBLE REPAIRS OR CHANGES ARE GUARANTEED IN THE ABOVE DESCRIBED WAY FOR THE REST OF THE ORIGINAL GUARANTEE TIME. THE GUARANTEE IS LIMITED TO THE REPAIR AND CHANGE OF THE FAULTY PRODUCT.

(THERE ARE NO OTHER IMPLICIT OR EXPLICIT GUARANTEES. MOREOVER M.B. WILL NOT BE RESPONSIBLE FOR CASUAL, CONSEQUENTIAL OR SPECIAL DAMAGES OR FOR OTHER DAMAGES, COSTS OR EXPENSES EXCEPT THE ONES FOR REPAIR OR CHANGE AS ABOVE DESCRIBED).

## VIII. DECLARATION OF CONFORMITY

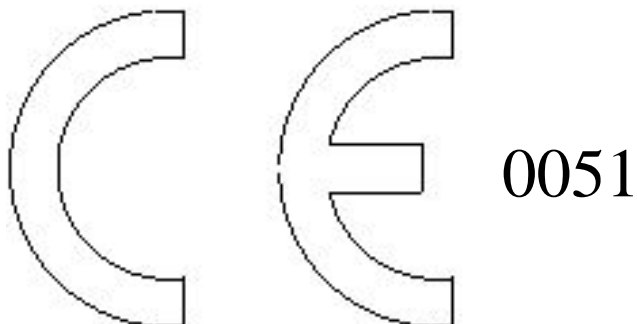Company M.B. PRO.MEC. SRL declares, under its responsibility, that the equipment called "DrySter", cod. MB-STERxx-AS000, is in compliance with requirement of emergency of D.P.R. 459/96, to directive EMC 89/336/CEE, directive 93/42/EEC, norm CEI EN 61010-1, norm CEI EN 61010-2-40 and all the norms and integrations to they connected.

1. Category: Dry Sterilizer
2. Company: M.B. PRO.MEC. SRL
3. Type: xx Lt
4. Series number:
5. Construction year: 2009
6. Date: 15/07/2009
7. Signature